


Sterilization tapes

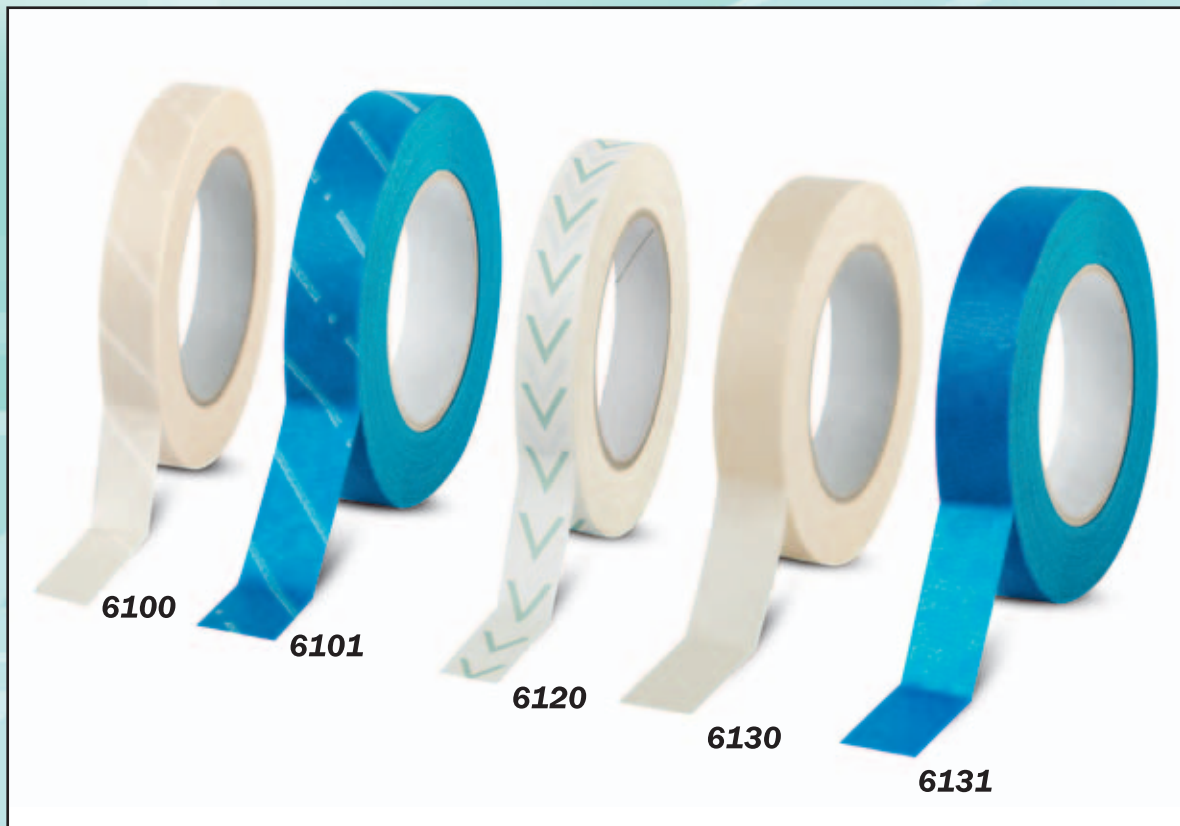

Sterilization tapes range is specially formulated for packing purposes as well as to show visibly that the tape has undergone a sterilization cycle.


The backing of the different tapes consists of strong semicreped impregnated paper which gives excellent conformability during application.

Sterilization tapes unwind easily and can be easily torn by hand.

Special adhesives formulated to withstand high temperature and humidity ensure excellent bond to paper, plastic, non-wovens, board, metal, glass and linen.


When removed after sterilization, specific tapes do not leave any glue residue.

The indicator inks on the surface of the tape change colour after exposure to its proper sterilization method: steam or dry heat.

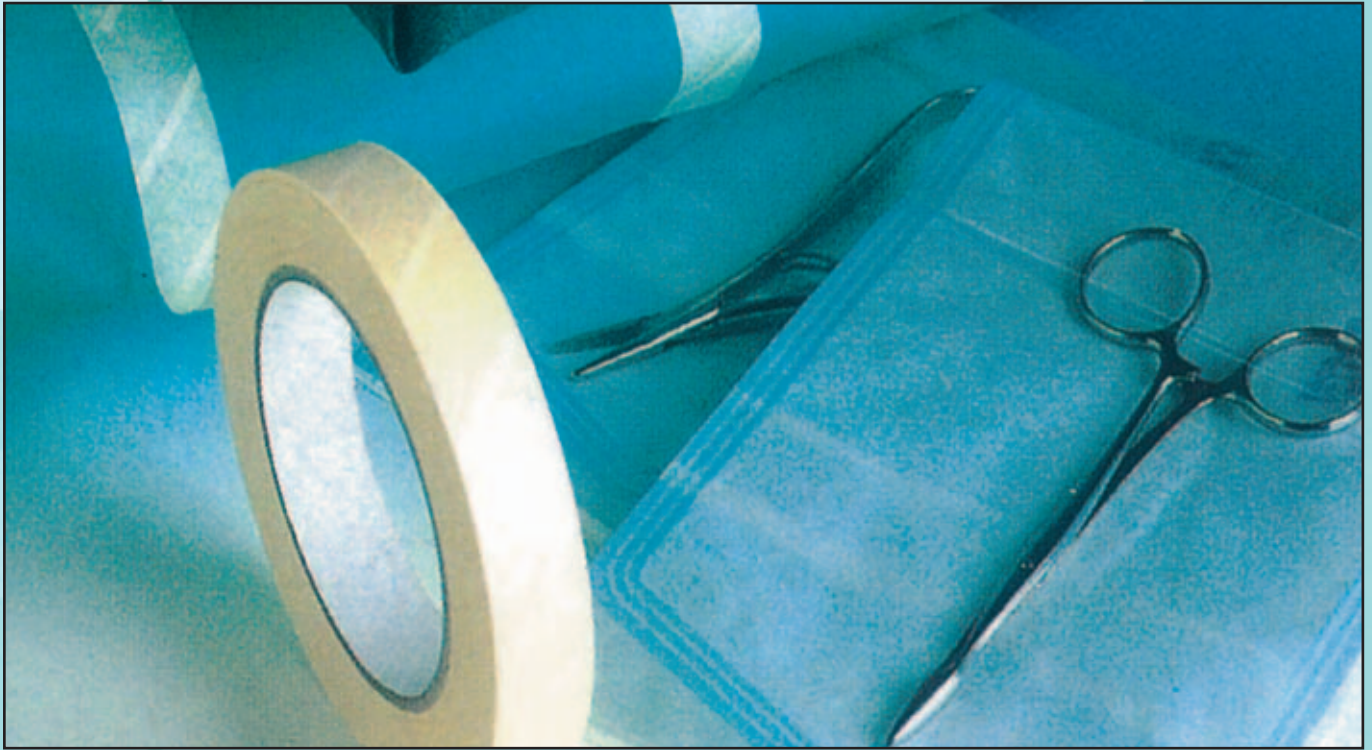

TECHNICAL TAPES RANGE FOR AUTOCLAVE STERILIZATION

	6100	6101	6120	6130	6131
BACKING	WHITE SATURATED PAPER	BLUE SATURATED PAPER	WHITE SATURATED PAPER	WHITE SATURATED PAPER	BLUE SATURATED PAPER
TOTAL THICKNESS	0,135 mm	0,160 mm	0,135 mm	0,135 mm	0,160 mm
ADHESION	4 N/25 mm	5,5 N/25 mm	8 N/25 mm	4 N/25 mm	5,5 N/25 mm
ELONGATION	7,0%	7,0%	7,5%	7,0%	7,0%
STERILIZATION TERMS	saturated steam cycles from 121°C to 134°C	saturated steam cycles from 121°C to 134°C	40' a 160°C	saturated steam cycles from 121°C to 134°C	saturated steam cycles from 121°C to 134°C
COLOUR INDICATOR	Before: White After: Brown	Before: White After: Brown	Before: Green After: Brown	None	None

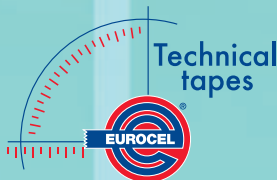

SICAD GROUP

Via Caduti Liberazione 57 - 21040 Uboldo (Va) ITALIA
Tel. 0039 02 964121 r.a. - Fax 0039 02 96780832 - 96412322

E-mail: info@eurocel.it - www.eurocel.it

SICAD DO BRASIL LTDA

AV. Pio XII N°1388 - B. Porto Alegre CEP 13360 - 000 - Capivari - San Paolo
Fone 0055/19/3492-9710 - Fax 0055/19/3491-2672

E-mail: sicad@sicad.com.br - www.eurocel.it

SICAD FRANCE

Rue de la Zamin - BP 232 Capinghem
59462 Lomme Cedex - France
Tél. 0033 03 20170017 - Fax 0033 03 20936652

E-mail: contact@sicadfrance.fr - www.eurocel.it

SICAD SOUTH AFRICA (PTY) LTD

South Africa P.O. Box 1500 - Germiston 1400
Tel. +27 (0)11 873-1311 - Fax +27 (0)11 873-2396

E-mail: rburns@sicad.co.za - www.eurocel.it



*IMPORTANT

The above mentioned data represent the averages obtained in our Laboratory. Suggestions for use and final results correspond to our best technical experience but cannot be intended as our commitment of liability. The user will determine the suitability of the method of applications for each specific use.